

4.0 Home Page Navigation

After sign in, the user is navigated to the homepage (*figure 4.1*) which provides access to the following features:

- 1) Cycle selection and execution (**Section 4.1**)
- 2) Cycle deletion (**Section 4.2**)
- 3) Manual sound controls (**Section 4.3**)
- 4) Manual light controls (**Section 4.4**)
- 5) Cycle creation (**Section 5**)
- 6) System Sign out (**Section 6**)

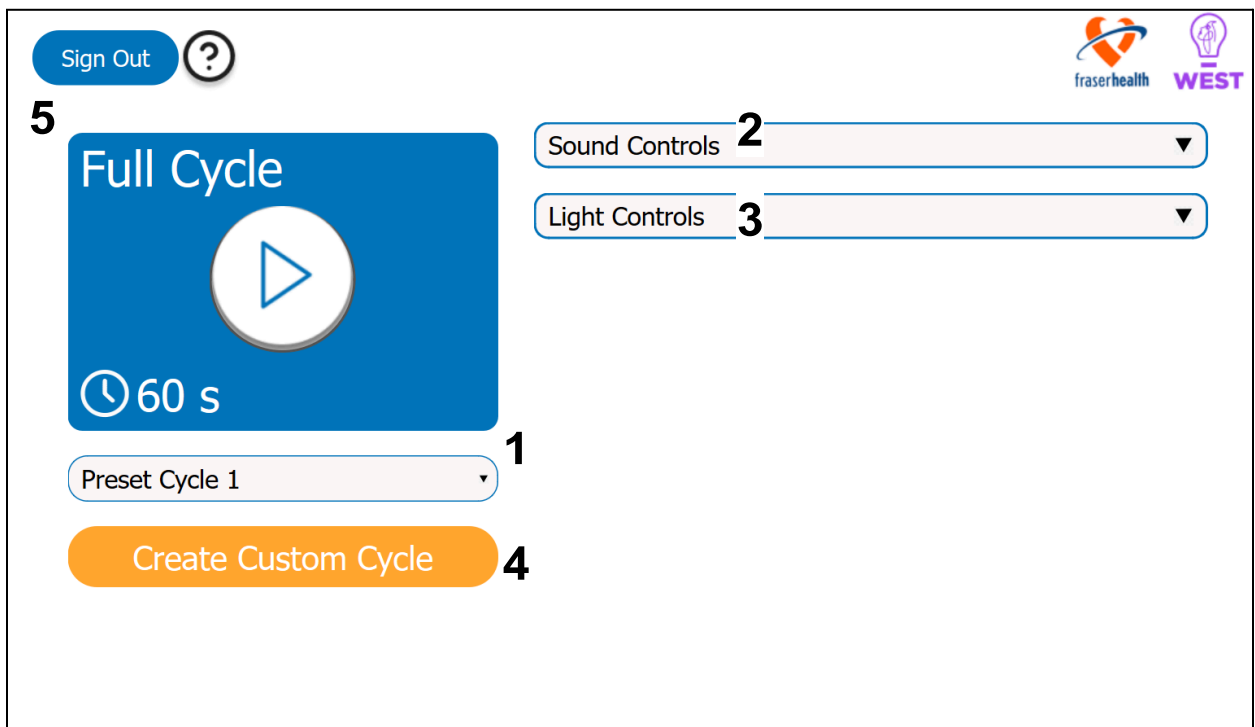


Fig. 4.1 Home Screen

4.1 Cycle Selection and Execution

A “Cycle” is a collection of sounds with a set brightness played over a set duration that when executed simulates the experience of an MRI machine. The system initially comes with 3 non-deleteable preset cycles for simulation, however custom cycles can be created and loaded as seen in **Section 5: Create Custom Cycle**. These cycles can be selected and played on the MRI simulator using the cycle selection drop down menu shown by *figure 4.1, pointer 1*.

In order to select a cycle the user can touch the dropdown menu to view the collections of cycles saved to the system. Pressing on the desired cycle will select it, preparing it to be executed. As a reference “Preset Cycle 2” is selected and is highlighted in *figure 4.2*.

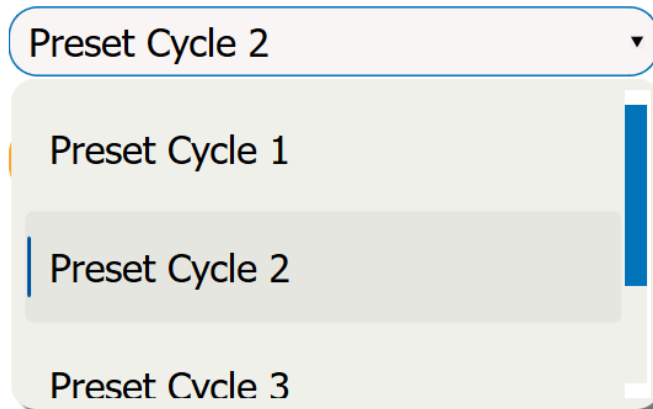


Fig. 4.2 Cycle Drop Down Menu

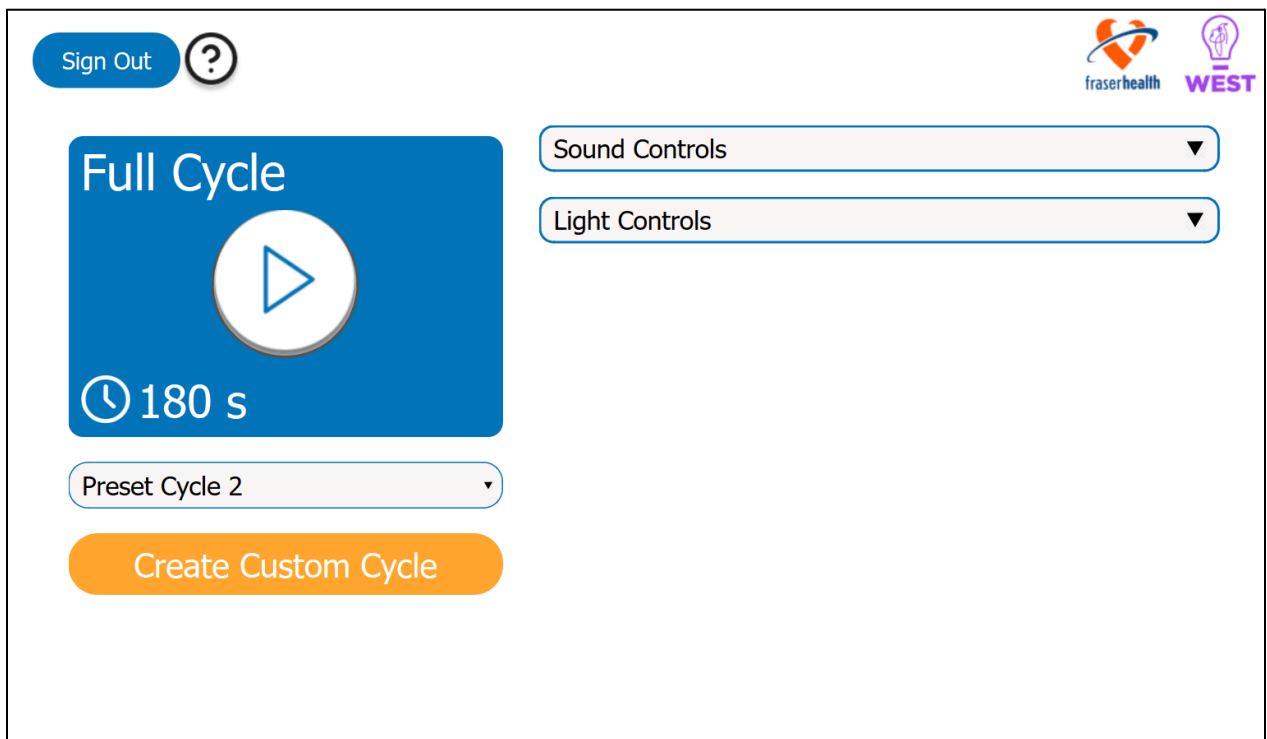


Fig. 4.3 Home Page With “Preset Cycle 2 Selected”

Once a cycle is selected it will be displayed on the home page as seen in *figure 4.3*. The large cycle play button, pictured top left in the blue square, can then be pressed to play the selected cycle. This will reveal a summary of the cycle contents (*figure 4.4*)

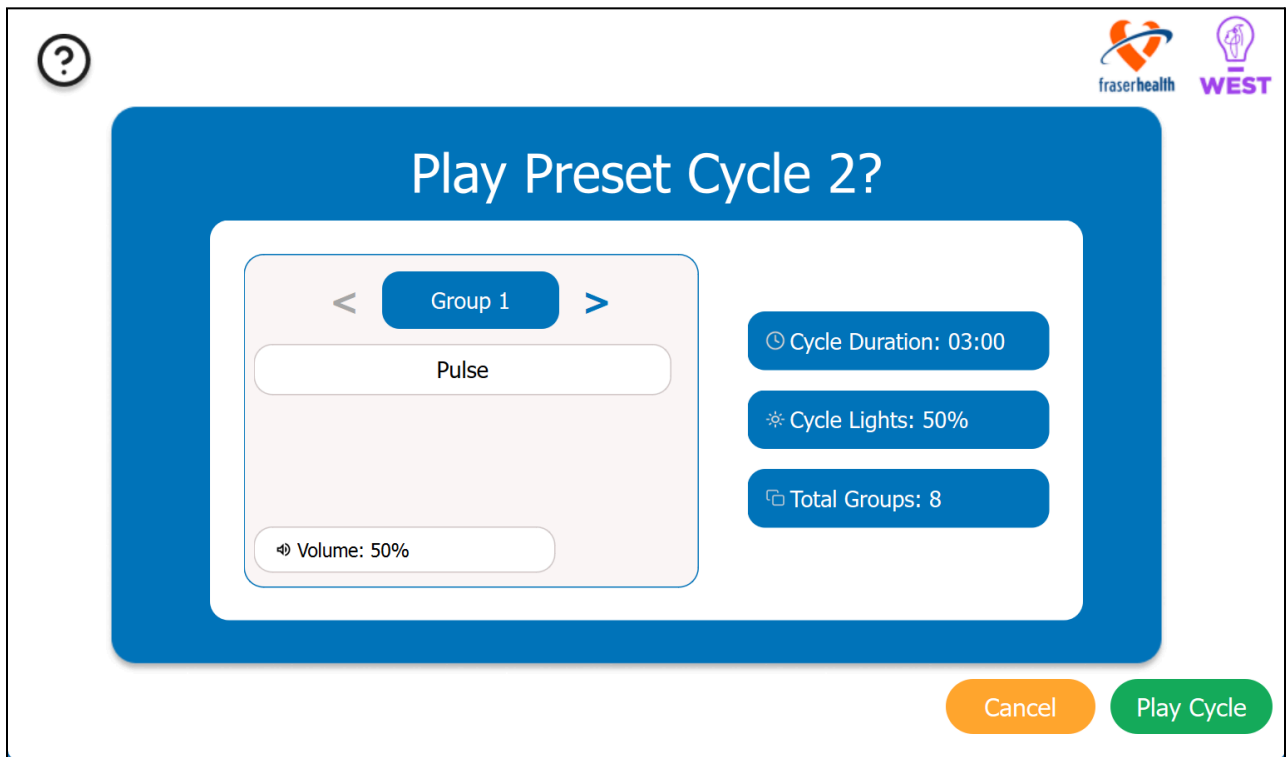


Fig. 4.4 Cycle Summary Page

The cycle summary page explains the contents of the selected cycle. Attributes such as duration, cycle brightness, and total sound groups are displayed on the right. In cycles, sounds are played concurrently and at the same volume per group. These groups play in succession, providing variation of sounds to simulate the experience of an MRI. Viewing the contents of these groups can be done using the “groups” panel, found on the left. Navigating across groups can be done using blue arrows near the top of the panel.

The user will now either cancel cycle execution using the “cancel” button, returning the user to the homepage, or continue to cycle execution by pressing the “play cycle” button in the bottom right. This will begin playing the cycle, revealing a countdown for the remaining duration of the cycle as seen in *figure 4/5*.

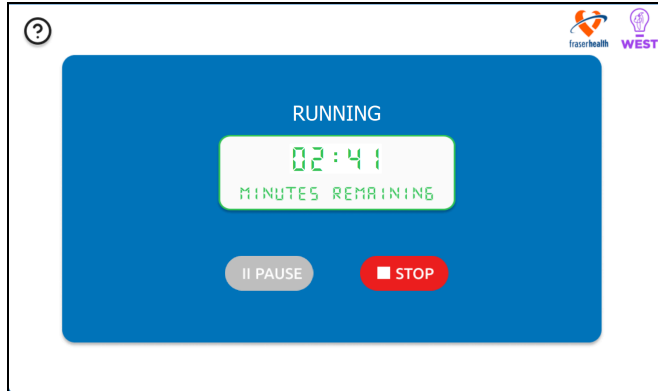


Fig. 4.5 Cycle Running Page

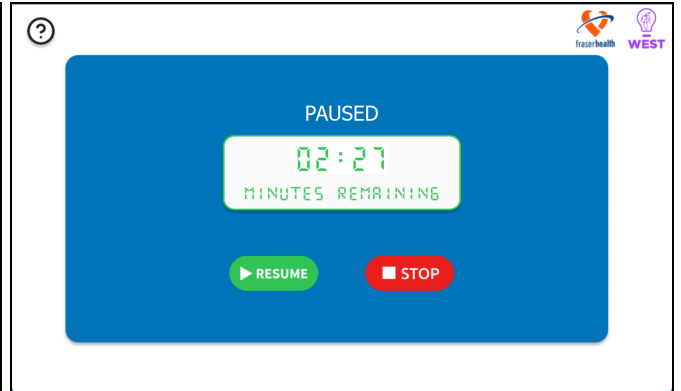


Fig. 4.6 Cycle Paused Page

As the cycle plays the user may press the “PAUSE” button, which pauses the cycle. This can then be resumed using the “RESUME” button later on as seen on *figure 4.6*. If the user presses the “STOP” button, the cycle execution immediately stops and returns the user to the home page. If the cycle completes in its entirety, the user is also automatically returned to the home page.

4.2 Cycle Deletion

Cycle deletion removes a cycle from system memory. Deletion can only be performed on cycles not labelled as “Preset”, namely cycle 1, 2, and 3. Using the example of cycle 4 in *figure 4.7*, after pressing the cycle drop down menu, any cycles which have the trash can  icon next to them may be deleted. Deleting a cycle will reveal a confirmation page, where the user can confirm or cancel deletion of a cycle.

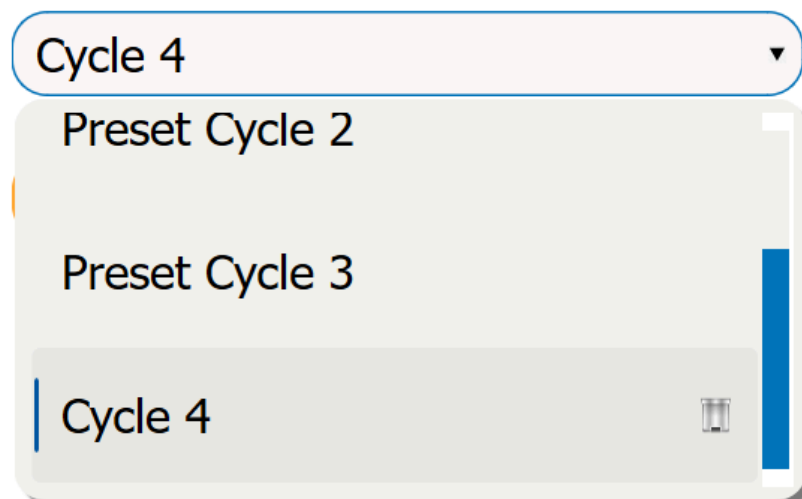


Fig. 4.7 Cycle Selection Dropdown

4.3 Manual Sound Controls

When the sound control panel (*figure 4.8*) is pressed, the panel expands and eight sounds will be available (*figure 4.9*). Tapping a sound once will play the sound aloud. Clicking the same sound again will stop the sound from playing. A total of 3 sounds can be played concurrently.

Directly below the sound selection is a slider adjustment for sound volume. The volume can be increased or decreased by clicking the “plus” and “minus” buttons respectively. Note that this panel will immediately collapse when accessing other features.



Fig. 4.8 Sound Control Panel (closed)

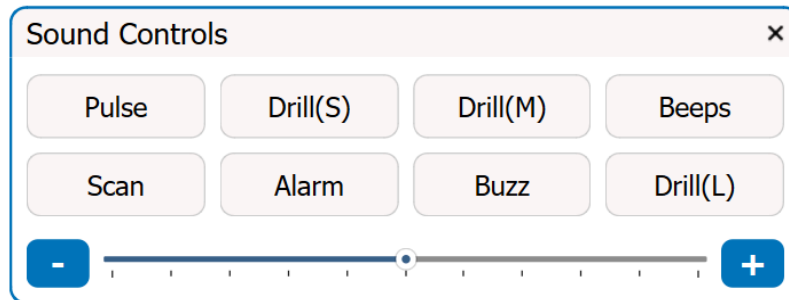


Fig. 4.9 Sound Control Panel (open)

4.4 Manual Light Controls

The brightness of the MRI simulator can be controlled by pressing the “Light Controls” panel. The default setting is set to medium intensity, and can be increased or decreased by clicking the “plus” and “minus” buttons respectively. They can also be toggled on and off using the “OFF/ON” button. Note that this panel will immediately collapse when accessing other features.

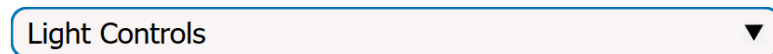


Fig. 4.10 Light Control Panel (closed)

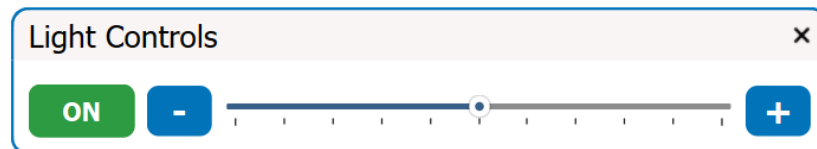


Fig. 4.11 Light Control Panel (open)